



**yuwell**



Technical Specifications

Model		
	YF-01	YN-02
Type	Full Face Mask	Nasal Mask
Weight	119g	64g
Therapy Pressure Range	4-30 cmH2O	4-30 cmH2O
Tube Connection	φ22mm	φ22mm
Flow Resistance (At 50 L/min)	≤0.5cmH2O	≤0.5cmH2O
Flow Resistance (At 100 L/min)	≤1.5cmH2O	≤1.5cmH2O
Expiratory Resistance with Anti Asphyxia Valve Open (At 50 L/min)	≤2.0cmH2O/(L/S)	—
Inspiratory Resistance with Anti Asphyxia Valve Open (At 50 L/min)	≤2.0cmH2O/(L/S)	—
Anti Asphyxia Valve Open Pressure	≤2.0cmH2O/(L/S)	—
Anti Asphyxia Valve Closed Pressure	≤2.0cmH2O/(L/S)	—

JIANGSU YUYUE MEDICAL EQUIPMENT & SUPPLY CO., LTD

Sales Center Address: Huanyuan East Road No.1, Xuzhuang Software Park, Nanjing, Jiangsu Province, P.R. China

Tel: 86-25-8558 2701 Fax: 86-25-5832 2068

Manufacturer: Suzhou Yuyue Medical Technology Co., Ltd.

Manufacturing Address: No.9 Jinfeng Road, Suzhou Science & Technology Town, 215163 Suzhou, Jiangsu, PEOPLE'S REPUBLIC OF CHINA

Tel: 86-512-67373001

Email:sales@yuyue.com.cn & yuwellsales@yuyue.com.cn

www.yuwell.com



YUWELL  
BREATHWEAR MASK SERIES  
CLASSIC MODELS

2020, Yuyue medical, all rights reserved. The "YUYUE", yuwell, the yuwell logo and name are trademarked by Jiangsu Yuyue Medical Equipment & Supply Co., Ltd, and may not be distributed, modified, reproduced in whole or in part without the prior. All specifications and product configurations are subject to change without notification.

www.yuwell.com

# YUWELL BREATHWEAR MASK SERIES

## CLASSIC MODELS



Fit



comfort



convenience



light & soft




### Rebound Silicone Cushion

The mask is made of rebound silicone cushion, which is airtight, light and soft. It is designed to reduce facial pressure, improve fit and comfort.




### Adaptive Silicone Cushion

The silicone cushion has an adaptive design at the nose bridge, which can be automatically adjusted according to different face shape, better relieve the pressure of the high nose and improve user experience.

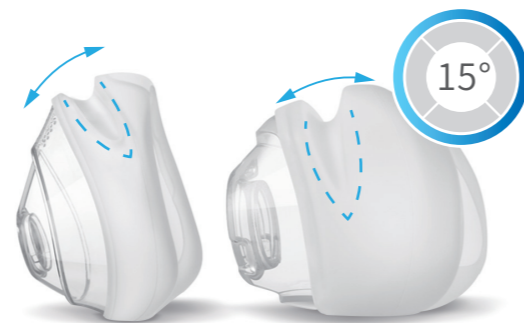

### Quiet Vent-hole

Honeycomb design, reducing airflow erosion intensity, gentle and quiet.

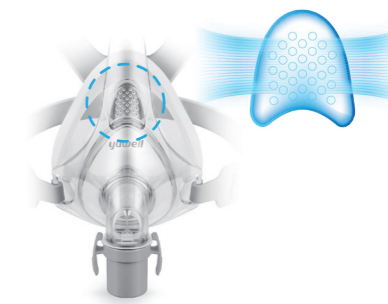

### Quick-release Connector

Using a quick release connector easily separate or connect the mask and the device.



### Soft And Light Headgear

Selected soft skin-friendly fabrics, light and stable, help patients adapt to treatment comfortably.

